



LIMS and workflow integration as impulse to establish Virtual Microscopy to daily clinical use

K. Saeger¹, S. Wienert¹, K. Polus², D. Emmel³, P. Hufnagl⁴

¹VMscope GmbH, Berlin

²Nexus/Dis GmbH, Frankfurt a.M.

³Loescap Technology GmbH, Berlin

⁴Charite Universitätsmedizin Berlin



Current situation of VM in clinical use


In Clinical applications, Virtual Microscopy is used for:

- Tele-consultation: Second opinion, supervision of remote pathologies
- Frozen-section telepathology (controversial)
- Case presentations: clinical conferences, continuing education
- Documentation and archiving of rare cases

...but not for the daily clinical pathology.



Advantages

- High degree on automation: Saves costs
- Digital archive: Saves space and time
- Remote work: Faster diagnosis for rare cases
- Support through image analysis:
Faster diagnosis and more reproducibility
- More diagnostic functionality:
Staining overlay, tele-consultation, ...

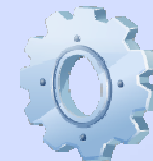


Disadvantages

- Additional process step: Digitalisation
- Huge hardware investments
- Training of personnel
- Diagnosis on the monitor is unfamiliar and legally indistinct

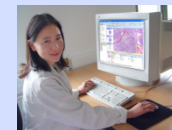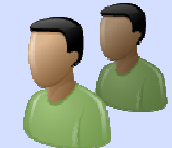


Already technically realised

- Unsupervised, automatic mass scanning
- Automatic slide registration through barcodes
- Visualisation of slides on the monitor
- Annotating virtual slides and creating snapshots

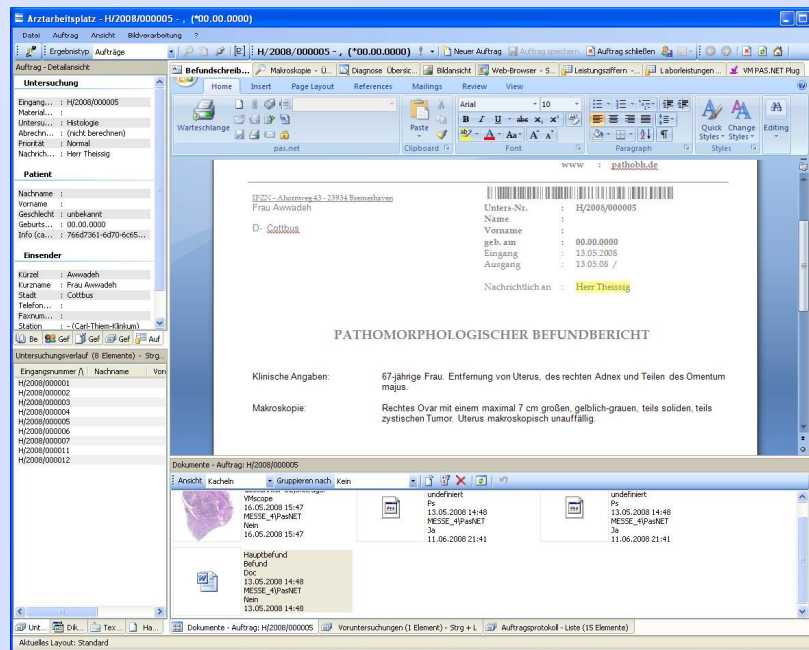



Current status of the digital workflow

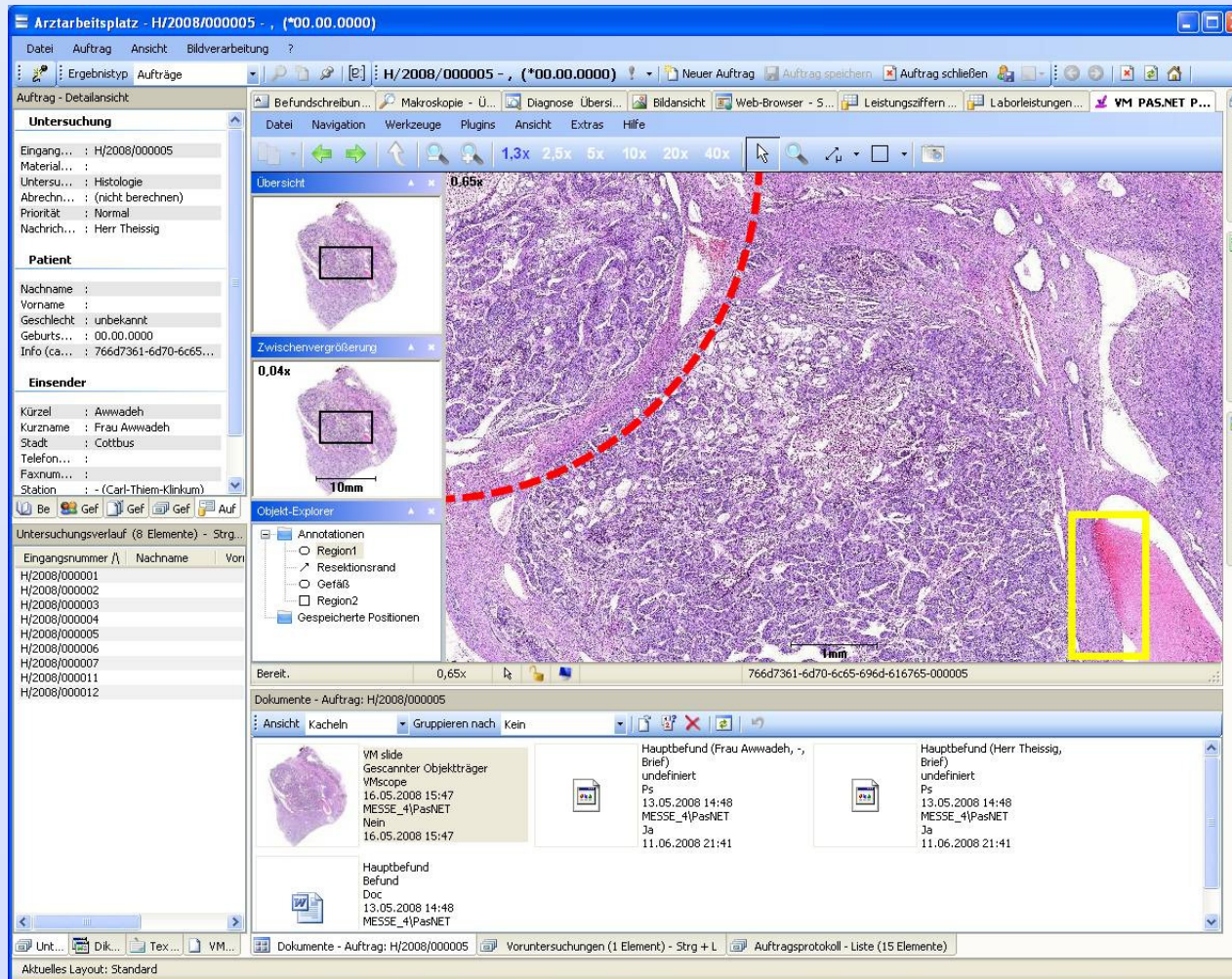

“VM Slide Explorer 2010” by VMscope GmbH



Remaining barriers for digital clinical pathology

Hardware automation

- Integrated autostainers with interface to LIMS/Pathology system
- Durability of glass slide barcodes during the laboratory processes
- Autostainers/coverslippers with hardware interface to slide scanners

Software integration

- LIMS integration
- PACS integration

Software automation and user support


- Workflow integration
- Specialised Applications (Tumour marker quantification)

A target LIMS/WSI solution would need:

- Digital laboratory order communication
- Slide scanner communication (preference of cito cases)
- System supervised auto-stainers and other laboratory machines
- Direct automated transfer of images to the storage in a DICOM-PACS
- Diagnosis on the virtual slide within the system
- Interfaces to the Virtual Microscope for:
 - Clinical reporting
 - Image analysis and
 - Tele-consultation systems


DICOM WG26/Supp.145 compatible interface to the clinic PACS with:

- Sharing of clinical short-term and long-term storage strategies
- Archiving of virtual slides of different slide scanners as WSI IOD
- WSI image streaming from the PACS
- Secured intra-clinic image communication
- Anonymous internet/tele-consultation image provision

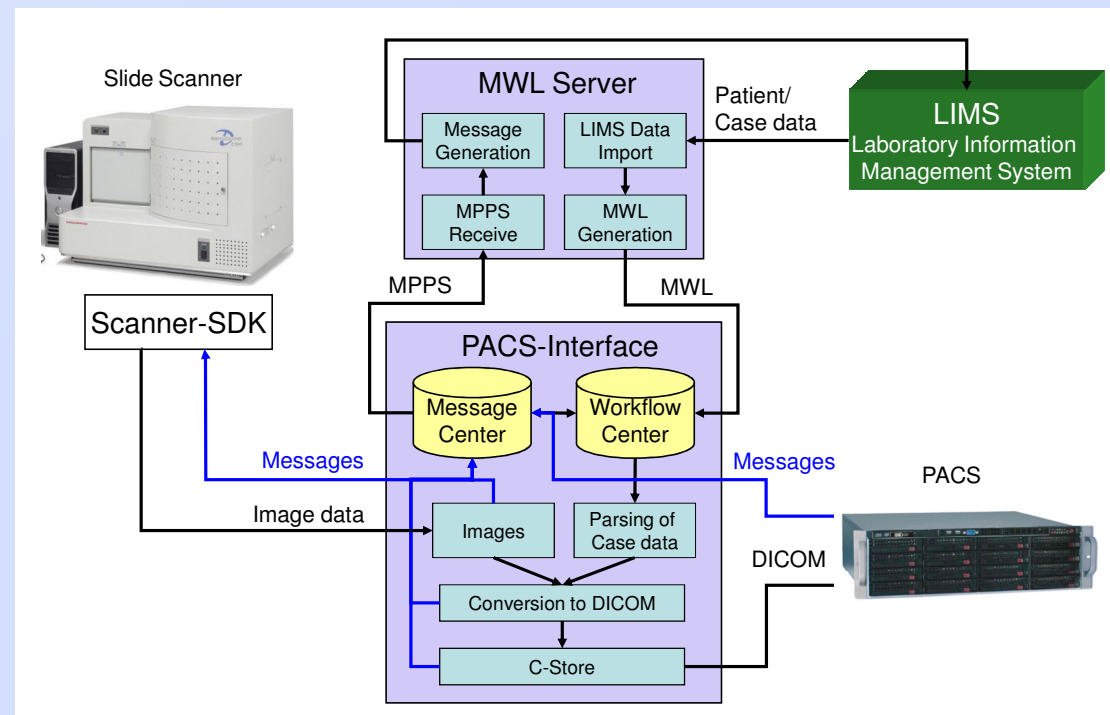

With a LIMS/PACS integration the workflow can be integrated and automated:


- Clinical pathways for decision for stainings, amount of slides, etc.
- Assignment to pathologists depending on case, duty schedule, and even personal workload.
- Organization of tele-consultation or remote case review by the LIMS

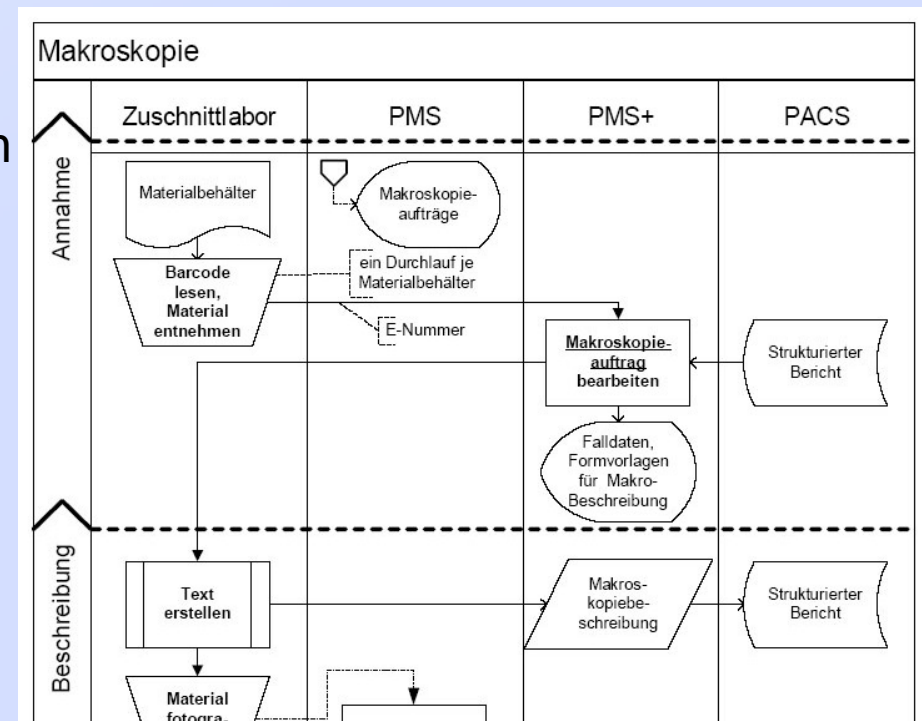



The user benefit increases with more specialised Applications:

- Tumour marker quantification
- Nucleus detection and evaluation

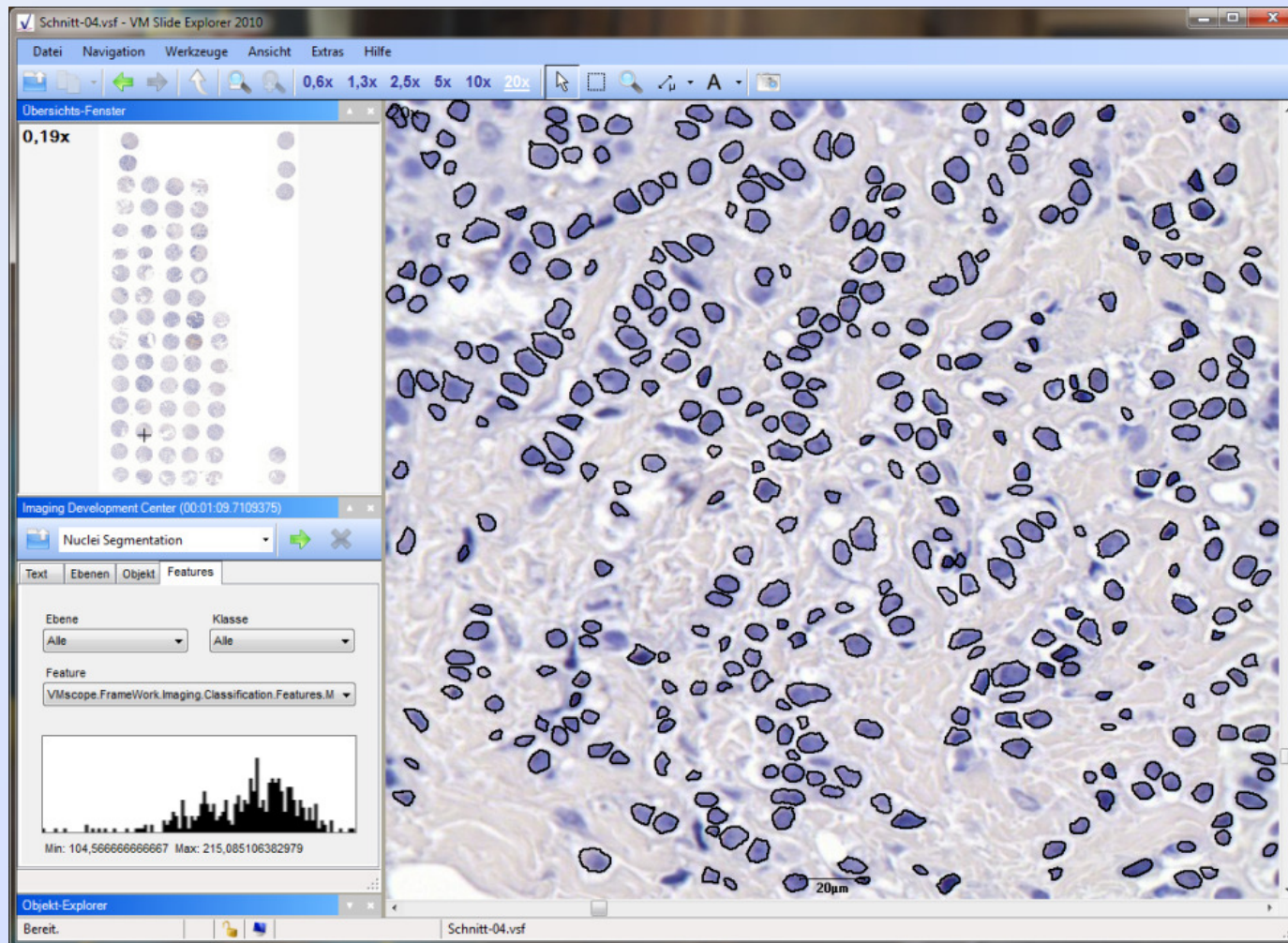

Automatic Nucleus Detection by VMscope GmbH



Thank you for your attention.

...and please visit our booth to learn more about VMscope.

This work has been supported by Zukunftsfonds Berlin (TSB, IBB).

